




With the unique VITA SYSTEM 3D-MASTER all natural tooth shades can be systematically determined and completely reproduced.

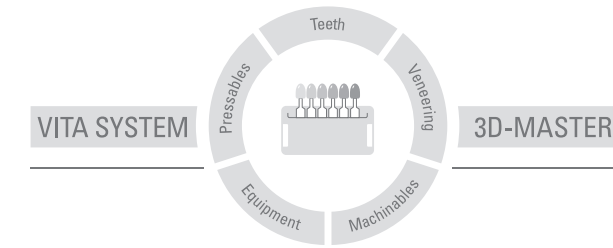

Please note: Our products must be used in accordance with the instructions for use. We accept no liability for any damage resulting from incorrect handling or usage. The user is furthermore obliged to check the product before use with regard to its suitability for the intended area of application. We cannot accept any liability if the product is used in conjunction with materials and equipment from other manufacturers that are not compatible or not authorized for use with our product. Furthermore, our liability for the accuracy of this information is independent of the legal basis and, in as far as legally permissible, shall always be limited to the value as invoiced of the goods supplied, excluding value-added tax. In particular, as far as legally permissible, we do not assume any liability for loss of earnings, indirect damages, ensuing damages or for third-party claims against the purchaser. Claims for damages based on fault liability (culpa in contrahendo, breach of contract, unlawful acts, etc.) can only be made in the case of intent or gross negligence. The VITA Modulbox is not necessarily a component of the product. Date of issue of this information: 09.14


After the publication of these information any previous versions become obsolete. The current version can be found at www.vita-zahnfabrik.com

VITA Zahnfabrik has been certified in accordance to the Medical Device Directive and the following products bear the CE
VITA V60 i-Line®

VITA

VITA Zahnfabrik H. Rauter GmbH & Co.KG
Spitalgasse 3 · D-79713 Bad Säckingen · Germany
Tel. +49(0)7761/562-0 · Fax +49(0)7761/562-299
Hotline: Tel. +49(0)7761/562-222 · Fax +49(0)7761/562-446
www.vita-zahnfabrik.com · info@vita-zahnfabrik.com
[facebook.com/vita.zahnfabrik](https://www.facebook.com/vita.zahnfabrik)

1989E - 0914 (0.5) Si - Version (02)

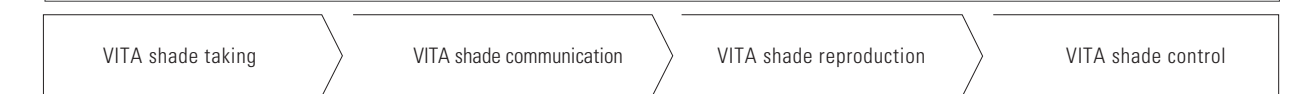

Date of issue 09.2014

VITA shade, VITA made. **VITA**

VITA V60 i-Line® – A true classic without compromises.



When it comes to dental ceramics and dental ceramics furnaces, VITA stands for top quality and proven technology. These characteristics are clearly evident in the VITA V60 i-Line. The VITA V60 i-Line focuses on two key aspects: ease of use and excellent firing results.

In developing this dental furnace, we decided to reduce the range of functions to the essentials and provide unbeatable value for the money. The VITA V60 i-Line has an impressive touch display, an intuitive user interface and the durable VITA fused quartz muffle for perfect firing results. The sturdy metal housing and the integrated finishing tables are extremely durable and resistant to damage.

The VITA V60 i-Line was developed on the basis of global customer surveys and decades of experience in the VITA research and development department. The result is a device that will meet all your expectations for a ceramics furnace.

VITA V60 i-Line® – Advantages



Excellent firing results



Very user-friendly basic, service and control programs

1. User-friendly: the touch display and intuitive menu navigation make working with the VITA V60 i-Line very easy. The minimalist design and robust technology focus on one central function: reliable firing.

2. Plug & play: no complicated assembly - just plug the furnace in, switch it on and start work!

3. Made in Germany: state-of-the-art engineering and top quality.

4. VITA quality: use of the best materials and components with just one goal: uncompromising reliability!

5. Easy to service: modular design, easy access to all wearing parts and well-organized arrangement make it easy to service and reduce your costs.

6. High performance ceramics furnace: a real workhorse. Extremely robust, integrated power surge protection and long-lasting VITA fused quartz muffle.

7. Reduced to the essentials: a classic and proven furnace design, which intentionally avoids costly additional functions. The design of this device focuses on an essential quality: absolute reliability.



VITA V60 i-Line® – Technical data and characteristics

Technical specifications

- Suitable for all brands of dental ceramics (including lithium disilicate ceramic and VITA SUPRINITY)
- Service and automatic control and monitoring programs such as automatic temperature calibration before the program starts (VITA AutoAdjust), automatic temperature and vacuum control as well as firing data memory
- Short-term compensation of power supply fluctuations
- Adjustable vacuum hold times
- Lift monitoring
- Individually adjustable lift positions for pre-drying and slow cooling processes
- Integrated rapid cooling program
- VITA Energy Efficiency for TÜV-certified energy efficiency

| Technical Data | |
|----------------------------|---|
| Width | 260 mm |
| Depth | 420 mm |
| Height | 570 mm |
| Weight | 18 kg |
| Housing | Sheet steel, powder coated |
| Firing chamber diameter | 90 mm |
| Firing chamber height | 55 mm |
| Maximum firing temperature | Max. 1200°C |
| Power supply | 230 Volt, 50 Hz or 100/110 volts, 50/60Hz |
| Maximum energy usage | Max. 1,5 kW |
| Classification | Safety Class 1 |

