


# VITA YZ<sup>®</sup> T – for high-strength substructures

Translucent. Highly accurate. Zirconia.

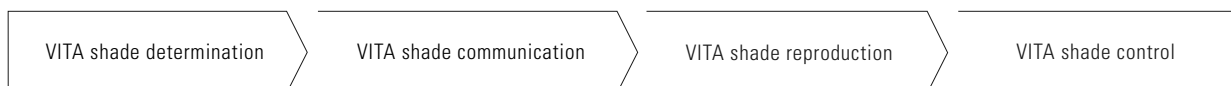

Date of issue: 10.16

VITA shade, VITA made.

**VITA**

## Product description

Thanks to its remarkable strength (approx. 1,200 MPa), VITA YZ T can be used to fabricate delicate substructures for crown and bridge restorations with a high precision of fit.

For exact shade reproduction, there is a system of matched components available for VITA YZ T, including coloring liquids, stains and veneering ceramic.

# VITA YZ<sup>®</sup> T – for high-strength substructures

## Indication

For substructures that offer high load capacity:

- crown substructures in the anterior/posterior areas
- multi-unit anterior and posterior bridge substructures

## Benefits

- high edge stability and precise sintering results for accurate restorations
- high reliability thanks to enormous strength and more than 10 years of clinical experience
- time-efficient fabrication thanks to high-speed sintering process of 80 minutes
- matched components (e.g. veneering ceramic) for exact shade reproduction

## Geometries

VITA YZ T blanks are available

- in block geometries: 20/19, 40/19, 55, 65/25
- in DISC geometries: Ø 98.4 in heights of 10, 12, 14, 16, 18, 20 and 25 mm

## Shades

VITA YZ T blanks are available

- as YZ T<sup>White</sup>: uncolored (white) for custom shading
- as YZ T<sup>Color\*</sup>: precolored in the shade/lightness levels Color LL1/light, LL2/medium und LL3/intense (only in heights of 14 and 18 mm)

\*) VITA YZ T<sup>Color</sup> block geometries are only available in LL1/light.

## System components



## VITA YZ® T blanks

### Block geometries, non-colored

Designation	Shade	Geometry	Contents	Prod. No.
VITA YZ T <sup>White</sup>	non-colored	YZ-20/19: 20 x 19 x 15.5 mm	4 pieces	EC4YW3007804
VITA YZ T <sup>White</sup>	non-colored	YZ-40/19: 39 x 19 x 15.5 mm	2 pieces	EC4YW3008042
VITA YZ T <sup>White</sup>	non-colored	YZ-55: 55 x 19 x 15.5 mm	2 pieces	EC4YW3008162
VITA YZ T <sup>White</sup>	non-colored	YZ-65/25: 65 x 25 x 22 mm	1 piece	EC4YW3008401

### DISC geometries, precolored

Designation	Shade	Geometry	Contents	Prod. No.
VITA YZ T <sup>Color</sup>	LL1/light	YZ-20/19: 20 x 19 x 15.5 mm	4 pieces	EC4YW3617804
VITA YZ T <sup>Color</sup>	LL1/light	YZ-40/19: 39 x 19 x 15.5 mm	2 pieces	EC4YW3618042
VITA YZ T <sup>Color</sup>	LL1/light	YZ-55: 55 x 19 x 15.5 mm	2 pieces	EC4YW3618162
VITA YZ T <sup>Color</sup>	LL1/light	YZ-65/25: 65 x 25 x 22 mm	1 piece	EC4YW3618401

### Disc geometries, non-colored

Designation	Shade	Geometry	Contents	Prod. No.
VITA YZ T <sup>White</sup>	non-colored	Ø 98.4 / h 10 mm	1 piece	ECDYW3981000
VITA YZ T <sup>White</sup>	non-colored	Ø 98.4 / h 12 mm	1 piece	ECDYW3981200
VITA YZ T <sup>White</sup>	non-colored	Ø 98.4 / h 14 mm	1 piece	ECDYW3981400
VITA YZ T <sup>White</sup>	non-colored	Ø 98.4 / h 16 mm	1 piece	ECDYW3981600
VITA YZ T <sup>White</sup>	non-colored	Ø 98.4 / h 18 mm	1 piece	ECDYW3981800
VITA YZ T <sup>White</sup>	non-colored	Ø 98.4 / h 20 mm	1 piece	ECDYW3982000
VITA YZ T <sup>White</sup>	non-colored	Ø 98.4 / h 25 mm	1 piece	ECDYW3982500

### DISC geometries, precolored

Designation	Shade	Geometry	Contents	Prod. No.
VITA YZ T <sup>Color</sup>	LL1/light	Ø 98.4 / h 14 mm	1 piece	ECDYC3981461
VITA YZ T <sup>Color</sup>	LL2/medium	Ø 98.4 / h 14 mm	1 piece	ECDYC3981462
VITA YZ T <sup>Color</sup>	LL3/intense	Ø 98.4 / h 14 mm	1 piece	ECDYC3981463
VITA YZ T <sup>Color</sup>	LL1/light	Ø 98.4 / h 18 mm	1 piece	ECDYC3981861
VITA YZ T <sup>Color</sup>	LL2/medium	Ø 98.4 / h 18 mm	1 piece	ECDYC3981862
VITA YZ T <sup>Color</sup>	LL3/intense	Ø 98.4 / h 18 mm	1 piece	ECDYC3981863
VITA YZ T <sup>Color</sup>	LL2/medium	Ø 98.4 / h 20 mm	1 piece	ECDYC3982062

## VITA YZ® T COLORING LIQUID

Set:

### VITA YZ® T COLORING LIQUID SET – Complete Set - Prod. No. EZ0CZ500

Quantity	Contents	Material
1 each	Bottle, 100 ml	VITA YZ T COLORING LIQUID light / pale, medium, intense, neutral
4	90 ml	Work containers
1	piece	Acrylic tweezers for VITA YZ T COLORING LIQUID
1	piece	VITA YZ T COLORING LIQUID shade indicator

Single-color sets:

**VITA YZ® T COLORING LIQUID SET light / pale** (prod. no. EZ0CZ67250)

Quantity	Contents	Material
1	Bottle 250 ml	VITA YZ T COLORING LIQUID light / pale
1	90 ml	Work container
1	piece	Acrylic tweezers for VITA YZ T COLORING LIQUID

**VITA YZ® T COLORING LIQUID SET medium** (prod. no. EZ0CZ68250)

Quantity	Contents	Material
1	Bottle 250 ml	VITA YZ T COLORING LIQUID medium
1	90 ml	Work container
1	piece	Acrylic tweezers for VITA YZ T COLORING LIQUID

**VITA YZ® T COLORING LIQUID SET intense** (prod. no. EZ0CZ69100)

Quantity	Contents	Material
1	Bottle 100 ml	VITA YZ T COLORING LIQUID intense
1	90 ml	Work container
1	piece	Acrylic tweezers for VITA YZ T COLORING LIQUID

**VITA YZ® T COLORING LIQUID SET neutral** (prod. no. EZ0CZ84100)

Quantity	Contents	Material
1	Bottle 100 ml	VITA YZ T COLORING LIQUID neutral
1	90 ml	Work container
1	piece	Acrylic tweezers for VITA YZ T COLORING LIQUID

**System compatibility**

**VITA YZ® T – SYSTEM SOLUTIONS**

VITA offers VITA YZ T with a specific holder system for the CAD/CAM system:

- inLab (Sirona Dental GmbH)

**VITA YZ® T – UNIVERSAL SOLUTIONS**

VITA is offering VITA YZ T in a universal disc geometry (Ø 98.4 mm) for the CAD/CAM systems:

- CORiTEC line (imes-icore GmbH)
- DMG ULTRASONIC line (DMG Mori AG)
- RXD line (Röders GmbH)
- inLab MC X5 (Sirona Dental GmbH)
- S1/S2/K3/K4/K5 Impression (vhf camfacture AG)

Note: VITA YZ T Discs can essentially be processed in all open CAD/CAM systems that can process discs with a diameter of 98.4 mm (including the circumferential groove).

**VITA**