

## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH

Revision date: 12.02.2020

Product code: 307

Page 1 of 9

#### SECTION 1: Identification of the substance/mixture and of the company/undertaking

##### 1.1. Product identifier

VITA ADIVA CERA ETCH

##### 1.2. Relevant identified uses of the substance or mixture and uses advised against

###### Use of the substance/mixture

Use as laboratory reagent

##### 1.3. Details of the supplier of the safety data sheet

Company name:	VITA Zahnfabrik H.Rauter GmbH & Co.KG	
Post-office box:	1338 79704 Bad Säckingen	
Telephone:	+49(0)7761-562-0	Telefax: +49(0)7761-562-299
e-mail:	info@vita-zahnfabrik.com	
Internet:	www.vita-zahnfabrik.com	

##### 1.4. Emergency telephone number:

+49-(0)761-19240

#### SECTION 2: Hazards identification

##### 2.1. Classification of the substance or mixture

###### Regulation (EC) No. 1272/2008

Hazard categories:  
 Substance or mixture corrosive to metals: Met. Corr. 1  
 Acute toxicity: Acute Tox. 2  
 Acute toxicity: Acute Tox. 3  
 Acute toxicity: Acute Tox. 4  
 Skin corrosion/irritation: Skin Corr. 1A  
 Serious eye damage/eye irritation: Eye Dam. 1  
 Hazard Statements:  
 May be corrosive to metals.  
 Fatal in contact with skin.  
 Toxic if swallowed.  
 Harmful if inhaled.  
 Causes severe skin burns and eye damage.  
 Causes serious eye damage.

##### 2.2. Label elements

###### Regulation (EC) No. 1272/2008

###### Hazard components for labelling

hydrofluorid acid  
 Sulphuric acid ... %

Signal word: Danger

###### Pictograms:



###### Hazard statements

H301	Toxic if swallowed.
H310	Fatal in contact with skin.
H314	Causes severe skin burns and eye damage.
H332	Harmful if inhaled.

## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH

Revision date: 12.02.2020

Product code: 307

Page 2 of 9

#### Precautionary statements

P405	Store locked up.
P260	Do not breathe dust/fume/gas/mist/vapours/spray.
P280	Wear protective gloves/protective clothing/eye protection/face protection.
P303+P361+P353	IF ON SKIN (or hair): Take off immediately all contaminated clothing. Rinse skin with water or shower.
P305+P351+P338	IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.
P310	Immediately call a POISON CENTER/doctor.

#### 2.3. Other hazards

No information available.

### SECTION 3: Composition/information on ingredients

#### 3.2. Mixtures

##### Hazardous components

CAS No	Chemical name			Quantity
	EC No	Index No	REACH No	
	GHS Classification			
7664-93-9	Sulphuric acid ... %			5 - < 10 %
	231-639-5	016-020-00-8		
	Skin Corr. 1A; H314			
7664-39-3	hydrofluorid acid			1 - < 5 %
	231-634-8	009-002-00-6	01-2119458860-33	
	Acute Tox. 1, Acute Tox. 2, Acute Tox. 2, Skin Corr. 1A; H310 H330 H300 H314			
64-17-5	ethyl alcohol			1 - < 5 %
	200-578-6	603-002-00-5	01-2119457610-43	
	Flam. Liq. 2, Eye Irrit. 2; H225 H319			

Full text of H and EUH statements: see section 16.

### SECTION 4: First aid measures

#### 4.1. Description of first aid measures

##### General information

First aider: Pay attention to self-protection! Remove affected person from the danger area and lay down.

##### After inhalation

Provide fresh air. Medical treatment necessary. IF INHALED: Call a POISON CENTER/doctor// if you feel unwell.

##### After contact with skin

After contact with skin, wash immediately with plenty of water and soap. Take off immediately all contaminated clothing and wash it before reuse. Call a physician immediately. IF ON SKIN (or hair): Ca-Gluconate solution Immediately call a POISON CENTER/doctor/.

##### After contact with eyes

In case of contact with eyes flush immediately with plenty of flowing water for 10 to 15 minutes holding eyelids apart and consult an ophthalmologist.Immediately call a POISON CENTER/doctor/.

##### After ingestion

Rinse mouth immediately and drink plenty of water. Do NOT induce vomiting. Adverse human health effects and symptoms: Gastric perforation. Call a physician immediately. If swallowed, immediately drink: Ca-Gluconate solution

## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH

Revision date: 12.02.2020

Product code: 307

Page 3 of 9

#### **4.2. Most important symptoms and effects, both acute and delayed**

No information available.

#### **4.3. Indication of any immediate medical attention and special treatment needed**

Treat symptomatically.

### SECTION 5: Firefighting measures

#### **5.1. Extinguishing media**

##### **Suitable extinguishing media**

Co-ordinate fire-fighting measures to the fire surroundings.

#### **5.2. Special hazards arising from the substance or mixture**

Non-flammable.

#### **5.3. Advice for firefighters**

Wear a self-contained breathing apparatus and chemical protective clothing. Full protection suit.

#### **Additional information**

Suppress gases/vapours/mists with water spray jet. Collect contaminated fire extinguishing water separately.

Do not allow entering drains or surface water.

### SECTION 6: Accidental release measures

#### **6.1. Personal precautions, protective equipment and emergency procedures**

Provide adequate ventilation. Do not breathe gas/fumes/vapour/spray. Avoid contact with skin, eyes and clothes. Use personal protection equipment.

#### **6.2. Environmental precautions**

Do not allow to enter into surface water or drains.

#### **6.3. Methods and material for containment and cleaning up**

Absorb with liquid-binding material (e.g. sand, diatomaceous earth, acid- or universal binding agents). Treat the recovered material as prescribed in the section on waste disposal.

#### **6.4. Reference to other sections**

Safe handling: see section 7

Personal protection equipment: see section 8

Disposal: see section 13

### SECTION 7: Handling and storage

#### **7.1. Precautions for safe handling**

##### **Advice on safe handling**

If handled uncovered, arrangements with local exhaust ventilation have to be used. Do not breathe gas/fumes/vapour/spray.

##### **Advice on protection against fire and explosion**

No special fire protection measures are necessary.

#### **7.2. Conditions for safe storage, including any incompatibilities**

##### **Requirements for storage rooms and vessels**

Keep container tightly closed. Keep locked up. Store in a place accessible by authorized persons only. Provide adequate ventilation as well as local exhaust at critical locations.

##### **Hints on joint storage**

No special measures are necessary.

#### **7.3. Specific end use(s)**

Use as laboratory reagent

### SECTION 8: Exposure controls/personal protection

## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH

Revision date: 12.02.2020

Product code: 307

Page 4 of 9

#### 8.1. Control parameters

##### Occupational exposure limit values

CAS No	Name of agent	ppm	mg/m <sup>3</sup>	fib/cm <sup>3</sup>	Category	Origin
7664-39-3	Hydrogen fluoride	1,8	1,5		TWA (8 h)	
		3	2,5		STEL (15 min)	
7664-93-9	Sulphuric acid (mist)	-	0,05		TWA (8 h)	

#### 8.2. Exposure controls



##### Appropriate engineering controls

If handled uncovered, arrangements with local exhaust ventilation have to be used. Do not breathe gas/fumes/vapour/spray.

##### Protective and hygiene measures

Remove contaminated, saturated clothing immediately. Draw up and observe skin protection programme. Wash hands and face before breaks and after work and take a shower if necessary. When using do not eat or drink.

##### Eye/face protection

Wear eye/face protection.

##### Hand protection

Recommended glove articles KCL Dermatril P Wearing time with occasional contact (splashes): 480 min  
When handling with chemical substances, protective gloves must be worn with the CE-label including the four control digits. The quality of the protective gloves resistant to chemicals must be chosen as a function of the specific working place concentration and quantity of hazardous substances. For special purposes, it is recommended to check the resistance to chemicals of the protective gloves mentioned above together with the supplier of these gloves.

##### Skin protection

Wear suitable protective clothing.

##### Respiratory protection

If handled uncovered, arrangements with local exhaust ventilation have to be used.

## SECTION 9: Physical and chemical properties

#### 9.1. Information on basic physical and chemical properties

Physical state:

Colour: light red  
Odour: characteristic

pH-Value: 2,0

##### Changes in the physical state

Melting point: not determined  
Initial boiling point and boiling range: 100 °C

##### Flammability

Solid: not applicable  
Gas: not applicable

Lower explosion limits: not determined

## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH

Revision date: 12.02.2020

Product code: 307

Page 5 of 9

Upper explosion limits:	not determined
<b>Auto-ignition temperature</b>	
Solid:	not applicable
Gas:	not applicable
Decomposition temperature:	not determined
<b>Oxidizing properties</b>	
Not oxidising.	
Vapour pressure: (at 50 °C)	<=1100 hPa
Density:	1,06000 g/cm <sup>3</sup>
<b>Solubility in other solvents</b>	
not determined	
Partition coefficient:	not determined
Vapour density:	not determined
Evaporation rate:	not determined

#### 9.2. Other information

Solid content:	1,5 %
----------------	-------

### SECTION 10: Stability and reactivity

#### 10.1. Reactivity

No hazardous reaction when handled and stored according to provisions.

#### 10.2. Chemical stability

The product is stable under storage at normal ambient temperatures.

#### 10.3. Possibility of hazardous reactions

No known hazardous reactions.

#### 10.4. Conditions to avoid

none

#### 10.5. Incompatible materials

No information available.

#### 10.6. Hazardous decomposition products

No known hazardous decomposition products.

### SECTION 11: Toxicological information

#### 11.1. Information on toxicological effects

##### **Acute toxicity**

Fatal in contact with skin.  
Toxic if swallowed.  
Harmful if inhaled.

##### **ATEmix calculated**

ATE (oral) 100,4 mg/kg; ATE (dermal) 100,4 mg/kg; ATE (inhalation vapour) 10,04 mg/l; ATE (inhalation aerosol) 1,004 mg/l

## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH

Revision date: 12.02.2020

Product code: 307

Page 6 of 9

CAS No	Chemical name				
	Exposure route	Dose	Species	Source	Method
7664-39-3	hydrofluorid acid				
	oral	ATE 5 mg/kg			
	dermal	ATE 5 mg/kg			
	inhalation vapour	ATE 0,5 mg/l			
	inhalation aerosol	ATE 0,05 mg/l			
64-17-5	ethyl alcohol				
	dermal	LD50 7060 mg/kg			

#### Irritation and corrosivity

Causes severe skin burns and eye damage.  
Causes serious eye damage.

#### Sensitising effects

Based on available data, the classification criteria are not met.

#### Carcinogenic/mutagenic/toxic effects for reproduction

Based on available data, the classification criteria are not met.

#### STOT-single exposure

Based on available data, the classification criteria are not met.

#### STOT-repeated exposure

Based on available data, the classification criteria are not met.

#### Aspiration hazard

Based on available data, the classification criteria are not met.

#### Additional information on tests

The mixture is classified as hazardous according to regulation (EC) No 1272/2008 [CLP]. Special hazards arising from the substance or mixture!

## SECTION 12: Ecological information

### 12.1. Toxicity

The product is not: Ecotoxic.

### 12.2. Persistence and degradability

The product has not been tested.

### 12.3. Bioaccumulative potential

The product has not been tested.

### 12.4. Mobility in soil

The product has not been tested.

### 12.5. Results of PBT and vPvB assessment

The product has not been tested.

### 12.6. Other adverse effects

No information available.

#### Further information

Do not allow to enter into surface water or drains. Do not allow to enter into soil/subsoil.

## SECTION 13: Disposal considerations

### 13.1. Waste treatment methods

#### Disposal recommendations

Do not allow to enter into surface water or drains. Do not allow to enter into soil/subsoil. Dispose of waste

## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH

Revision date: 12.02.2020

Product code: 307

Page 7 of 9



according to applicable legislation.

#### Contaminated packaging



This material and its container must be disposed of as hazardous waste. Handle contaminated packages in the same way as the substance itself.

### SECTION 14: Transport information



#### Land transport (ADR/RID)

<b>14.1. UN number:</b>	UN 2922
<b>14.2. UN proper shipping name:</b>	CORROSIVE LIQUID, TOXIC, N.O.S. (Wasser)
<b>14.3. Transport hazard class(es):</b>	8
<b>14.4. Packing group:</b>	II
Hazard label:	8+6.1
	 
Classification code:	CT1
Special Provisions:	274
Limited quantity:	1 L
Excepted quantity:	E2
Transport category:	2
Hazard No:	86
Tunnel restriction code:	E

#### Inland waterways transport (ADN)

<b>14.1. UN number:</b>	UN 2922
<b>14.2. UN proper shipping name:</b>	CORROSIVE LIQUID, TOXIC, N.O.S. (Wasser)
<b>14.3. Transport hazard class(es):</b>	8
<b>14.4. Packing group:</b>	II
Hazard label:	8+6.1
	 
Classification code:	CT1
Special Provisions:	274 802
Limited quantity:	1 L
Excepted quantity:	E2

#### Marine transport (IMDG)

<b>14.1. UN number:</b>	UN 2922
<b>14.2. UN proper shipping name:</b>	CORROSIVE LIQUID, TOXIC, N.O.S. (Wasser)
<b>14.3. Transport hazard class(es):</b>	8
<b>14.4. Packing group:</b>	II
Hazard label:	8+6.1
	 
Special Provisions:	274
Limited quantity:	1 L
Excepted quantity:	E2
EmS:	F-A, S-B

## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH


Revision date: 12.02.2020

Product code: 307

Page 8 of 9

#### Air transport (ICAO-TI/IATA-DGR)

<b>14.1. UN number:</b>	UN 2922
<b>14.2. UN proper shipping name:</b>	CORROSIVE LIQUID, TOXIC, N.O.S. (Wasser)
<b>14.3. Transport hazard class(es):</b>	8
<b>14.4. Packing group:</b>	II
Hazard label:	8+6.1

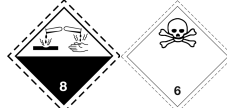

Special Provisions:	A3 A803
Limited quantity Passenger:	0.5 L
Passenger LQ:	Y840
Excepted quantity:	E2
IATA-packing instructions - Passenger:	851
IATA-max. quantity - Passenger:	1 L
IATA-packing instructions - Cargo:	855
IATA-max. quantity - Cargo:	30 L

#### 14.5. Environmental hazards

ENVIRONMENTALLY HAZARDOUS: no

#### 14.6. Special precautions for user

Warning: Toxic. strongly corrosive.

#### 14.7. Transport in bulk according to Annex II of Marpol and the IBC Code

not applicable

### SECTION 15: Regulatory information

#### 15.1. Safety, health and environmental regulations/legislation specific for the substance or mixture

##### National regulatory information

Employment restrictions:	Observe restrictions to employment for juveniles according to the 'juvenile work protection guideline' (94/33/EC). Observe employment restrictions under the Maternity Protection Directive (92/85/EEC) for expectant or nursing mothers.
Water hazard class (D):	2 - obviously hazardous to water
Skin resorption/Sensitization:	Permeates easily through outer skin and causes poisoning.

#### 15.2. Chemical safety assessment

Chemical safety assessments for substances in this mixture were not carried out.

### SECTION 16: Other information

#### Abbreviations and acronyms

ADR: Accord européen sur le transport des marchandises dangereuses par Route (European Agreement concerning the International Carriage of Dangerous Goods by Road)  
 IMDG: International Maritime Code for Dangerous Goods  
 IATA: International Air Transport Association  
 GHS: Globally Harmonized System of Classification and Labelling of Chemicals  
 EINECS: European Inventory of Existing Commercial Chemical Substances  
 ELINCS: European List of Notified Chemical Substances  
 CAS: Chemical Abstracts Service  
 LC50: Lethal concentration, 50%  
 LD50: Lethal dose, 50%



## Safety Data Sheet

according to Regulation (EC) No 1907/2006

### VITA ADIVA CERA ETCH

Revision date: 12.02.2020

Product code: 307

Page 9 of 9

#### Classification for mixtures and used evaluation method according to regulation (EC) No. 1272/2008 [CLP]

Classification	Classification procedure
Met. Corr. 1; H290	
Acute Tox. 2; H310	Calculation method
Acute Tox. 3; H301	Calculation method
Acute Tox. 4; H332	Calculation method
Skin Corr. 1A; H314	Calculation method
Eye Dam. 1; H318	Calculation method

#### Relevant H and EUH statements (number and full text)

H225	Highly flammable liquid and vapour.
H290	May be corrosive to metals.
H300	Fatal if swallowed.
H301	Toxic if swallowed.
H310	Fatal in contact with skin.
H314	Causes severe skin burns and eye damage.
H318	Causes serious eye damage.
H319	Causes serious eye irritation.
H330	Fatal if inhaled.
H332	Harmful if inhaled.

#### Further Information

The above information describes exclusively the safety requirements of the product and is based on our present-day knowledge. The information is intended to give you advice about the safe handling of the product named in this safety data sheet, for storage, processing, transport and disposal. The information cannot be transferred to other products. In the case of mixing the product with other products or in the case of processing, the information on this safety data sheet is not necessarily valid for the new made-up material.

*(The data for the hazardous ingredients were taken respectively from the last version of the sub-contractor's safety data sheet.)*